


**DIATOME**

**ultra sonic**

Thinner sections  
No compression  
Better structural preservation

Patented

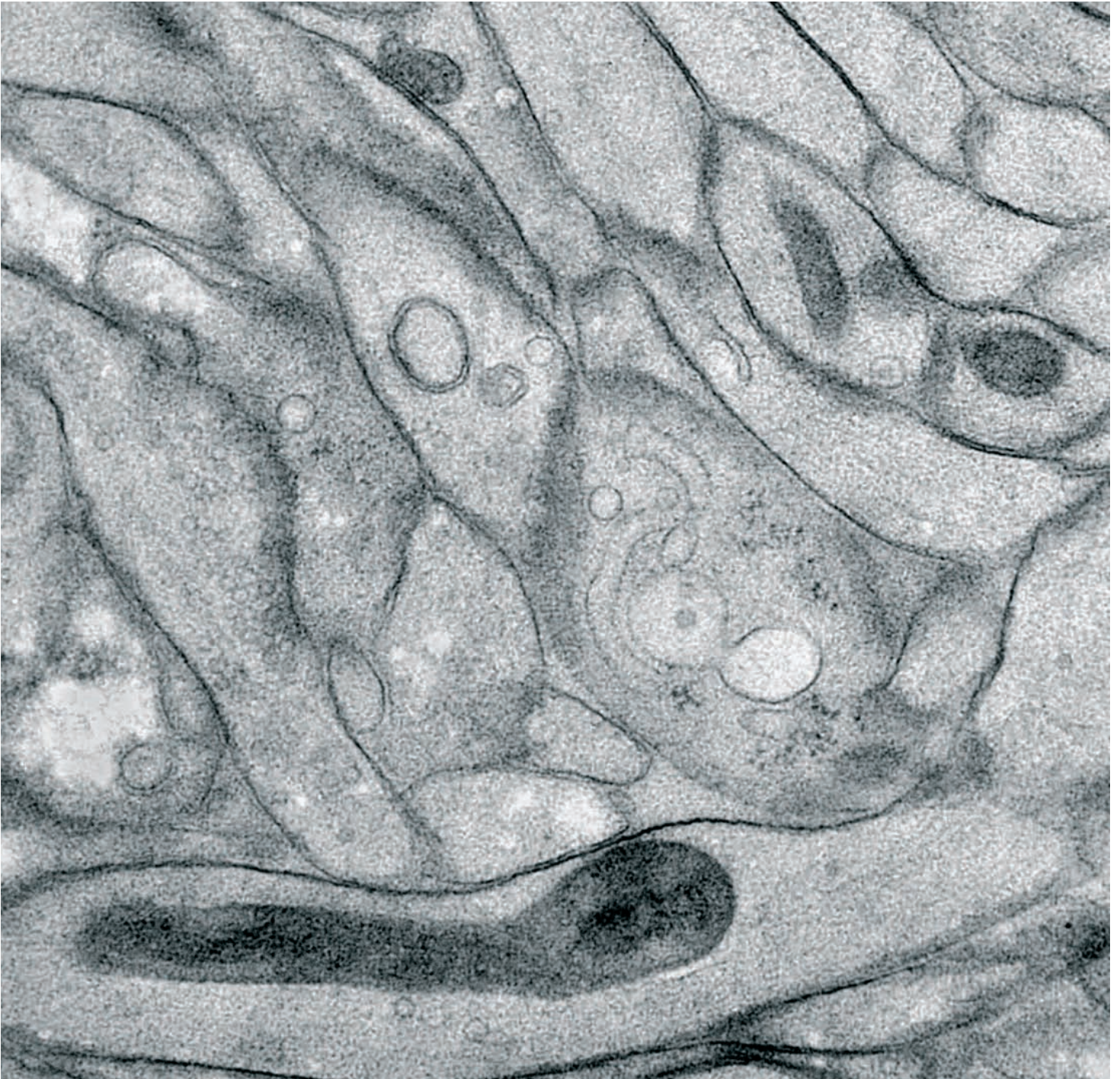





## The oscillating diamond knife, exclusively from Diatome!

For many years Diatome diamond knives have been used successfully for a wide range of room- and low temperature sectioning applications.

Major advances in immunocytochemistry, the sectioning of frozen hydrated specimens, semithin sectioning for optical microscopy, as well as the sectioning of hard industrial samples have been realised using Diatome diamond knives (refs. 1, 2, 3, 4, 5).

The development of the Static Line II ionizer enabled dry ultramicrotomy of Lowicryls (ref. 6) and a considerable improvement in cryosectioning (ref. 2).

Despite these innovations, until now, one major obstacle remained, preventing us from achieving perfect ultrathin sections: «compression» (we define compression as the shortening of the section compared to the sample height).

The amount of compression depends on various factors including:

- The wedge angle of the knife
- The hardness of the sample
- The interaction diamond surface / section surface (when sectioning dry)
- The section thickness.

The most critical factor is the wedge angle of the knife. It was shown that reducing the wedge angle results in a reduction of compression, hence better preservation of ultra-structure which allows a higher achievable resolution (refs. 7, 8).


However, the wedge angle may not be reduced ad infinitum. A further reduction results in a lower cutting edge quality and a considerably shorter life span.

In room temperature ultramicrotomy, we have found the following compression factors (section thickness 50nm):

- 10-20% for Epon, Araldite, EM-Bed, and other epoxy resins
- 12-24% for Lowicryl K4M
- 10-17% for Spurr's (hard grade)
- 8-13% for LR White (hard grade).

These limitations have stimulated our efforts to develop the oscillating knife.


Here we proudly introduce the Diatome ultra sonic, the oscillating diamond knife for room temperature ultramicrotomy. It was developed in collaboration with Dr Daniel Studer, Lab. of Anatomy, University of Berne.

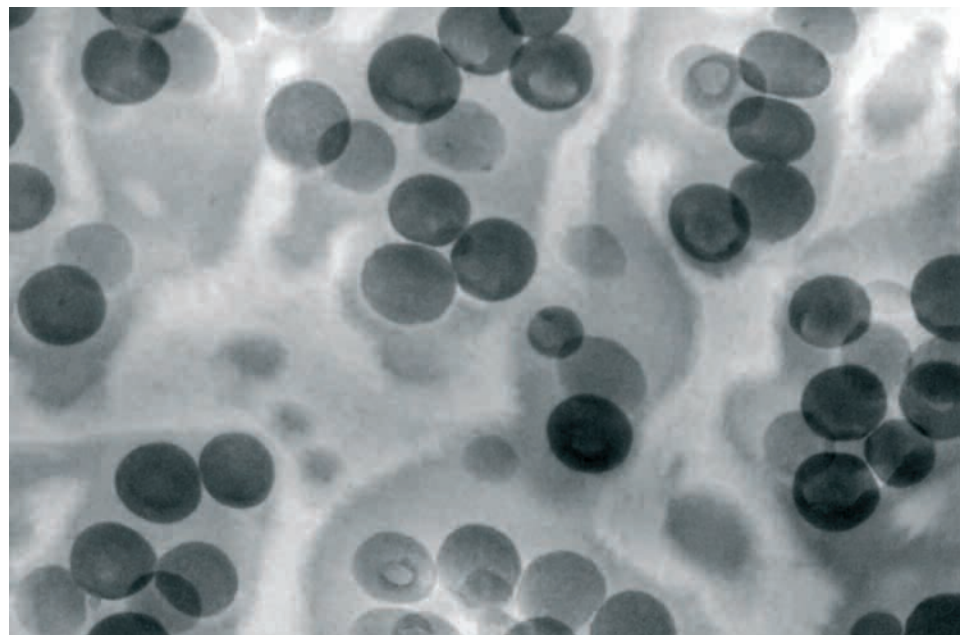

◀◀  
Neuropil of a HP frozen / freeze substituted hippocampus slice culture  
Daniel Studer and Werner Graber, Inst. of Anatomy, University of Bern  
└──────────┘ 200 nm

▲  
Polycarbonate modified with rubber  
Jens Sicking, Bayer Technology Services, Leverkusen

## ultra sonic

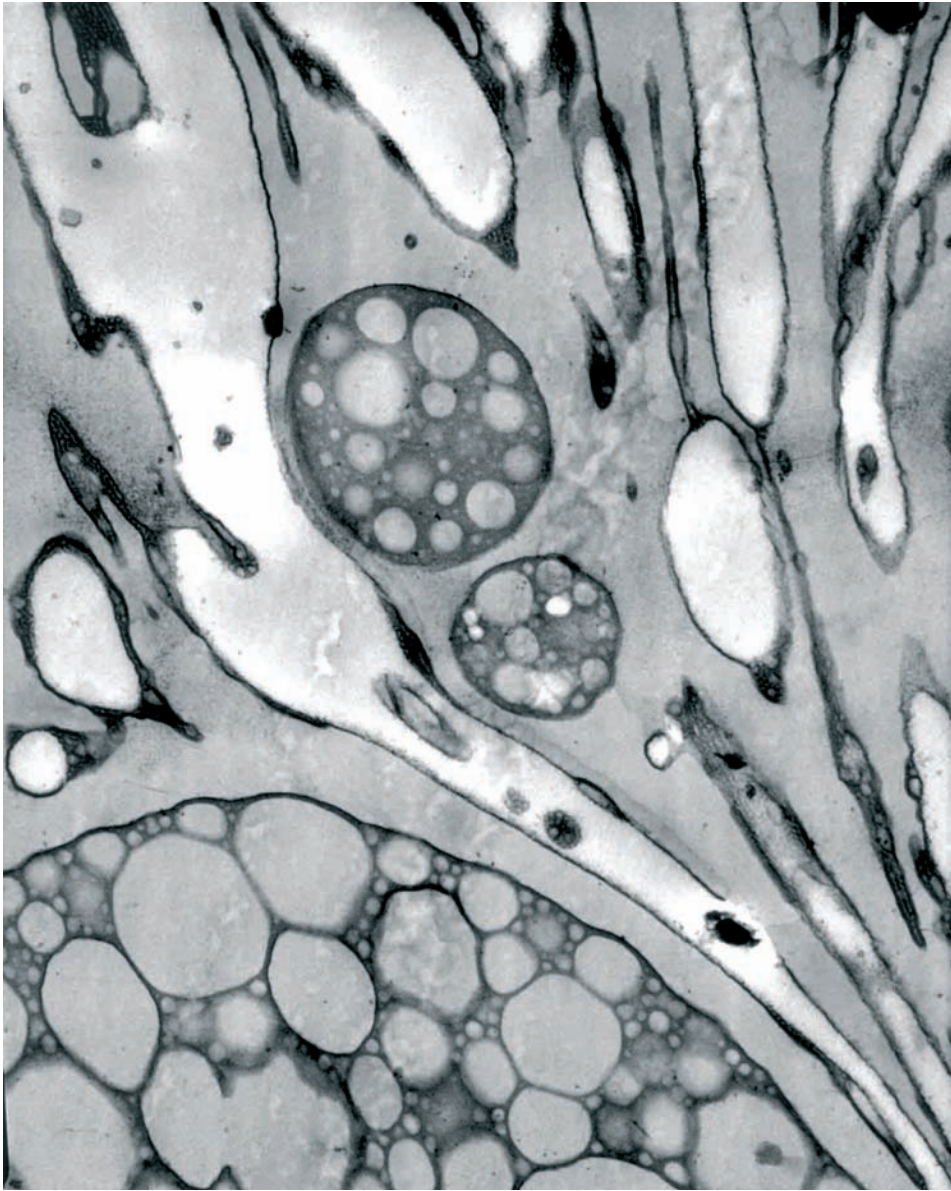

▲  
High Impact Polystyrene / Polypropylene blend  
Bob Vastenhou, Dow Chemicals, Terneuzen



A piezo actuator produces an oscillation of the knife at a desired frequency and amplitude, parallel with the cutting edge. A depression in the foot of the knife allows the oscillation parallel to the cutting edge. The depression is rigid in the north/south direction and guarantees stability in the cutting direction.


The ultra sonic knife produces ultrathin sections almost free of compression. The sections become thinner at the same thickness setting: since the volume of the section remains the same, the increased length leads to a decrease in thickness.

The theoretical considerations and first results of this invention were presented in the Journal of Microscopy (ref. 9).

Results with polymers using the oscillating knife were shown at the M&M 2002 meeting in Quebec (ref. 10) and in Microscopy Today (ref. 11).

We have tested the oscillating knife with the following samples:

- Biological samples in Epon, Araldite, EM Bed, etc.
- Biological samples in acrylic resins (Lowicryls, LR White)
- Rigid polymers such as PS, PMMA, ABS, HIPS, modified PP, etc.

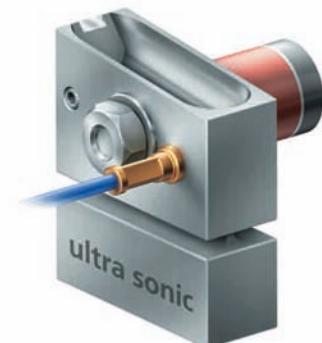

## Specifications

### Knife

Knife angle:	35°
Cutting range:	10-100nm
Cutting edge length:	3.0mm

### Control unit

Frequency:	20-45kHz, or automatic setting of the resonance
Amplitude:	variable (Voltage 0-30V)
Mains voltage:	230V, 110V

## References

1. H. Sitte: Advanced Instrumentation and Methodology related to Cryoultramicrotomy: a Review. Scanning Microscopy Supplement 10, pp. 87-466, 1996.
2. M. Michel, H. Gnägi and M. Müller: Diamonds are a cryosectioner's best friend. Journal of Microscopy, Vol. 166, Pt 1, pp. 43-56, 1992.
3. O.L. Reymond: The diamond knife „semi“: a substitute for glass or conventional diamond knives in the ultramicrotomy of thin and semi-thin sections. Bas. Applied Histochemistry, No. 30, pp. 487-494, 1986.
4. M.J.F. Blumer, P. Gahleitner, T. Narzt, C. Handl, B. Ruthensteiner: Strips of semi-thin sections: an advanced method with a new type of diamond knife. Journal of Neuroscience Methods 120, pp. 11-16, 2002.
5. Various publications on materials ultramicrotomy. Microscopy Research and Technique, Vol. 31. Number 4, pp. 265-310, 1995.
6. L. Edelmann: Freeze-substitution and the preservation of diffusable ions. Journal of Microscopy, Vol. 161, pp. 217-228, 1991.
7. J.C. Jésior: How to avoid compression. Journal of Ultrastructure and Molecular Structure Research, pp. 210-217, 1986.
8. J.C. Jésior: Use of low-angle diamond knives leads to improved ultrastructural preservation of ultrathin sections. Scanning Microscopy Supplement 3, pp. 147-153, 1989.
9. D. Studer and H. Gnägi: Minimal compression of ultrathin sections with use of an oscillating diamond knife. Journal of Microscopy, Vol. 197, Pt 1, pp. 94-100, 2000.
10. J.S. Vastenhout and H. Gnägi: Ultramicrotomy of polymers using an oscillating knife; improving polymer morphology. Microscopy and Microanalysis. 8 (Suppl. 2) 2002.
11. J.S. Vastenhout and J.D. Harris: improved sectioning of polymers using an oscillating diamond knife for Transmission Electron Microscopy. Microscopy Today, pp. 20-21, 2006.

Diatome US  
1560 Industry Road  
Hatfield PA 19440  
Phone 215 412 8390  
Fax 215 412 8450  
e-mail: sgkcck@aol.com  
www.emsdiasum.com